pangao

NECK MASSAGER

Instruction Manual

PG-2601B15

Shenzhen Pangao Electronic Co., Ltd.

Main Site: No.25 1st Industry Zone, Fenghuang Road,
Xikeng Village, Henggang Town, Longgang District,
Shenzhen, Guangdong China

Tel:+86-0755-33825988

Fax: +86-0755-33825989

Web: www.pan-go.com

Thank you for using Pangao product.

Please read this manual carefully before use and keep it well for future reference.
This manual includes technology and using information

Product Instruction

The PG-2601B15 intelligent neck massager is of a streamline ring-shaped design consistent with the human cervical vertebra physiological curvature. It combines four functions of far infrared heating, vibration, red light therapy and circular traction, which can lead the curative effect to reach directly the body's deep tissues and marrow. During the entire treatment period, it can always keep human body at a proper physical position to provide a scientific platform for correct prevention and treatment of cervical vertebra diseases, thus offering you an easy, efficient, cost-effective and pleasant multifunctional systemic health care at home.

Product Scope of Application

This massager can promote local blood circulation, relax local muscles and relieve fatigue, thus it's applicable to health care and prevention of neck pain, numbness of shoulders and arms, headache and dizziness arising from cervical vertebra diseases.

Product Features

- Ergonomic design, aesthetic and scientific, perfectly fit the neck curve;
- Synchronous health care system combining FIR heat, vibration, red light and circular traction is applied to form an efficient composite energy field;
- The FIR heat, which is designed with three grades "Low, Medium and High" levels, users may choose according to their needs and ambient temperature;
- It can promote local blood circulation with visible red light illumination with a wavelength of 600-650nm.
- Built-in high capacity lithium battery, USB fast charging, easy to use.
- LED display, easy to operate.
- Auto turn off without any operation after 30 minutes.

Main structure of product

USB cable

How to use

1. Please full charge the product before using it for the first time.
2. During the use of the massager, when you hear the buzzer beep for 3 consecutive times, it means that the built-in lithium battery power is low. You can plug one end of the charging cable into the host and plug another end into the mobile power bank or computer USB port (5V DC) charging, and the power indicator will flash green. When the battery is fully charged, the power indicator keeps green.

How to Operate


1. Press “” button shortly and the massager starts to work with default heating mode 2 and it will stop working automatically after 30 minutes.
2. Then press the “” button again to select 4 different modes, and the “1, 2, 3,4” indicators will light up accordingly. Cycle in sequence: 2 (medium heat)-3 (intermittent vibration + medium heat) - 4 (high heat) -1 (low heat) - 2(medium heat) - 3 (intermittent vibration + medium heat) - 4 (high heat) - 1 (low heat)...
3. Long press “” button 1~2 seconds and shut down the power.

Use of rechargeable battery

1. During the use of the massager, when you hear the buzzer beep for 3 consecutive times, it means that the built-in lithium battery power is low. You can plug one end of the charging cable into the host and plug another end into the mobile power bank or computer USB port (5V DC) charging, and then the power indicator will flash green. When the battery is fully charged, the power indicator keeps green. It will take about 2-3 hours to fully charge.
2. In the case of normal use, the longevity of this rechargeable battery is over 500 times for charging and discharging, and can be used for about 2-3 years (30 minutes per day).
3. When the rechargeable battery is not used for a long time, the service life will be shortened. It is recommended that it should be charged at least once a month.

The way to take out built-in rechargeable battery

1. The lithium battery built into this product contains substances that may pollute the environment. Please take out the battery before disposing of this product.
2. Please make sure to turn off the power before disassembling the product and taking out the lithium battery.
3. Follow the below steps from 1 to 4 to take out the battery gradually as the picture shown.
4. The battery must be discarded at a place regulated by the government or a recycling unit.

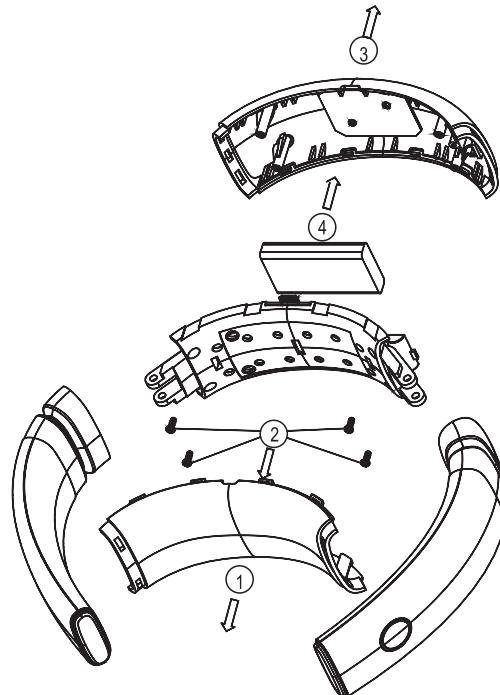Safety Precautions

1. Please read and understand all the contents of the manual before using this product.
2. This apparatus cannot be used for pregnant women and those in menstrual period and people with sensitive skin, heart disease, abnormal blood pressure, malignant tumor, cerebral vascular disease, or acute disease, or others under treatment before consultation with the doctor.
3. Prohibit using it when bathing, sweating or sleeping.
4. Don't use this massager on children or those without expressing ability consciousness.
5. Prohibit using it for patients with perceptual disability or those who are not sensitive to heat.
6. In case of any discomfort in using this massager, stop using immediately and seek a doctor.
7. Prohibit using it in high temperature, electromagnetic radiation, inflammable or damp places.
8. Please don't disassemble, repair and refit this massager without permission, or a trouble or electric shock will be caused.
9. Please don't wash this product with water.

Cleaning and Maintenance

1. When cleaning the product, wipe it with a soft cloth or a towel dampened with a little water, and do not wipe with a volatile preparation.
2. Please store it in a dry and ventilated place.
3. Don't store it in places with direct sunshine, high temperature or moisture.
4. Avoid using corrosive liquids such as benzene, thinner and other chemical solvents on this product and accessories.

Product scrapping

In order to avoid environmental pollution, when the massage device and accessories are aging and cannot be used, please do not discard them at random. Dispose of them. Dispose of the treatment according to local environmental protection requirements.

Troubleshooting

Do not disassemble or attempt repairs privately.

Problems	Inspection Item	Solution
Failure to boot	Is the internal rechargeable battery of the main unit low?	If it is in low power, please charge it with a USB charging cable.
Can't charge or no charging indication	Is the USB cable connected well?	Check that the USB cable is in place.
	Is the charging power supply low?	Is the charging power output normal?
Not working	in standby?	Press the model button again.

Technical Parameters and Specifications

1. Product specifications

Name: Neck massager
Model: PG-2601B15
Power supply: 5V DC 1A(USB);
Input power: 5W
Lithium battery: 3.7V DC 1100mAh
Heating temperature: +38°C~+48°C;
Charging time 2~3 hours
Dimension of packaging box: 210x188x77 mm
Weight: about 155g

2. Transport and storage conditions

Ambient temperature: -20°C ~ +55°C;
Relative humidity: ≤ 93%;
Atmosphere pressure: 50 kPa ~ 106 kPa.

NEKMASSAGEAPPARAAT

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PG-2601B15

Shenzhen Pangao Electronic Co., Ltd.

Hoofdlocatie: No.25 1st Industry Zone, Fenghuang
Road, Xikeng Village, Henggang Town, Longgang
District, Shenzhen, Guangdong China

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