## (EN) Water flosser



**Please note:** Read the instructions before use. Clean the device before using it for the first time by rinsing it thoroughly.

### Instructions:

1. Remove the reservoir from the pump and fill the reservoir with water or mouthwash up to the line. Now carefully put back the reservoir on the pump. **Always keep the water flosser upright to prevent water from leaking from the reservoir.**
2. Press the button a couple of times until water emerges from the mouth. You can adjust the speed of the flow by pressing the button hard or softly.
3. In order to clean your teeth, braces, implants, etc., hold the mouth against your teeth, on the edge of your gums. Press the button and move the flosser along your gums and between your teeth. Close your lips slightly to prevent splatters, but allow the water to run from your mouth into the basin freely.

### Maintenance:

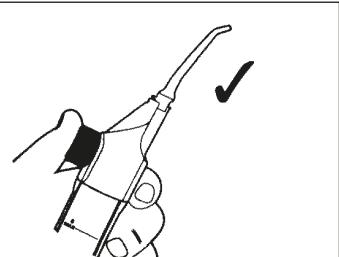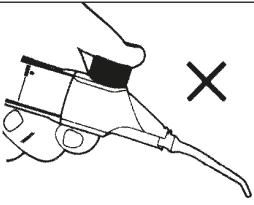1. Empty, clean and dry the reservoir after every use. Press the button to remove any remaining water from the pump.
2. Do not allow others to use this product for reasons of hygiene.
3. Keep in a safe place where it cannot get damaged.

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